


PhysioPedal[®]

User Guide

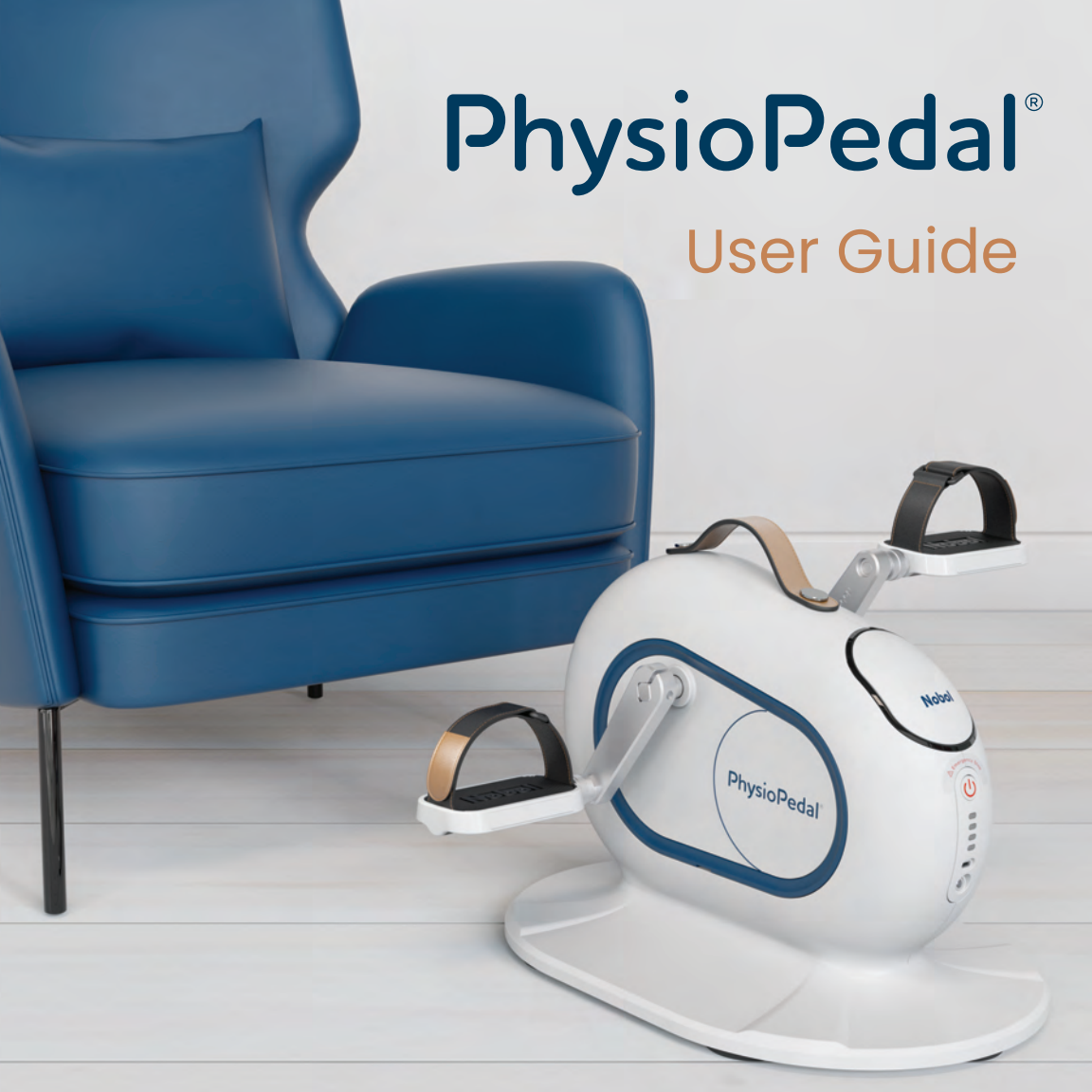

by **Nobel**[®]

**Congratulations
on your new
PhysioPedal**

CONTENTS

What's in the box	5
Safety instructions	6-7
PhysioPedal	8
Specs	9
Features	10-11
The display	12
The remote	13
Getting started	14-17
How to use	18-21
Cleaning and Storage	22
Warranty	23

WHAT'S IN THE BOX



- ① Body
- ② Battery
- ③ Wireless Remote
- ④ Charging Cable

- ⑤ Base
- ⑥ Pedals x 2
- ⑦ Handles x 2
- ⑧ Base Screws x 4

IMPORTANT SAFETY INSTRUCTIONS

When using The PhysioPedal, basic precautions should always be followed, including the following:

Read all instructions before using The PhysioPedal.

WARNING

To Reduce the risk of burns, fire, electric shock, or injury to persons:

Disconnect the battery pack from the cycle when not in use, and before putting on or taking off parts.

Use The PhysioPedal only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.

SAFETY INSTRUCTIONS

- Never operate The PhysioPedal if it's damaged, or not working properly, if it has been dropped or damaged, or dropped into water.
- Return The PhysioPedal to a service center for examination and repair.
- Never drop or insert any object into any opening.
- Do not use it in the rain.
- Do not operate where aerosol (spray) products are used or where oxygen is administered.
- For household use only.

SAVE THIS BOOKLET

PhysioPedal



LED Screen Display



Remote Control



Portable Battery



Patent Pending
Quick Release Pedal

SPECS



Power: 72W (Max)

Display: Digital

Battery Capacity: 2400mAh


Weight: 7.5Kg | 16Lbs

Battery Life: up to 3 hours

Output: 26V=1A

Input: 110~220V AC

Speed: 20~70 RPM


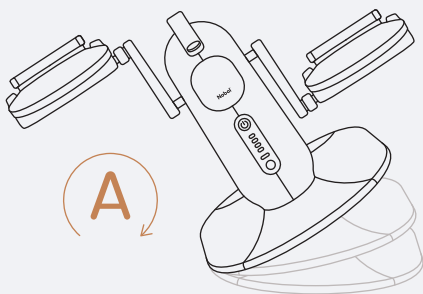


Built in safety

If the unit tilts it automatically stops.


Battery

High capacity rechargeable lithium-Ion battery.

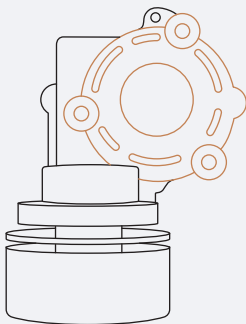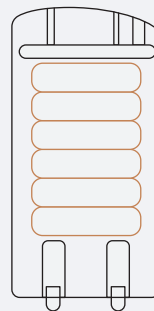

Motor

Medical grade brushless DC electric motor.


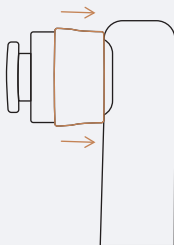


Quick Release Pedal

Pull the shaft and easelly insert or release the pedal.

Interactive Screen


Full tracking and reporting of all usage, calories, distance and time used.

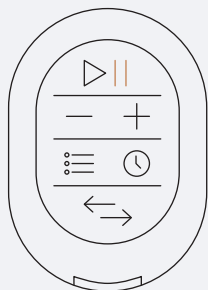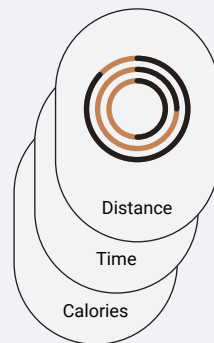


Remote Control

Senior friendly Intuitive Interface.


DISPLAY

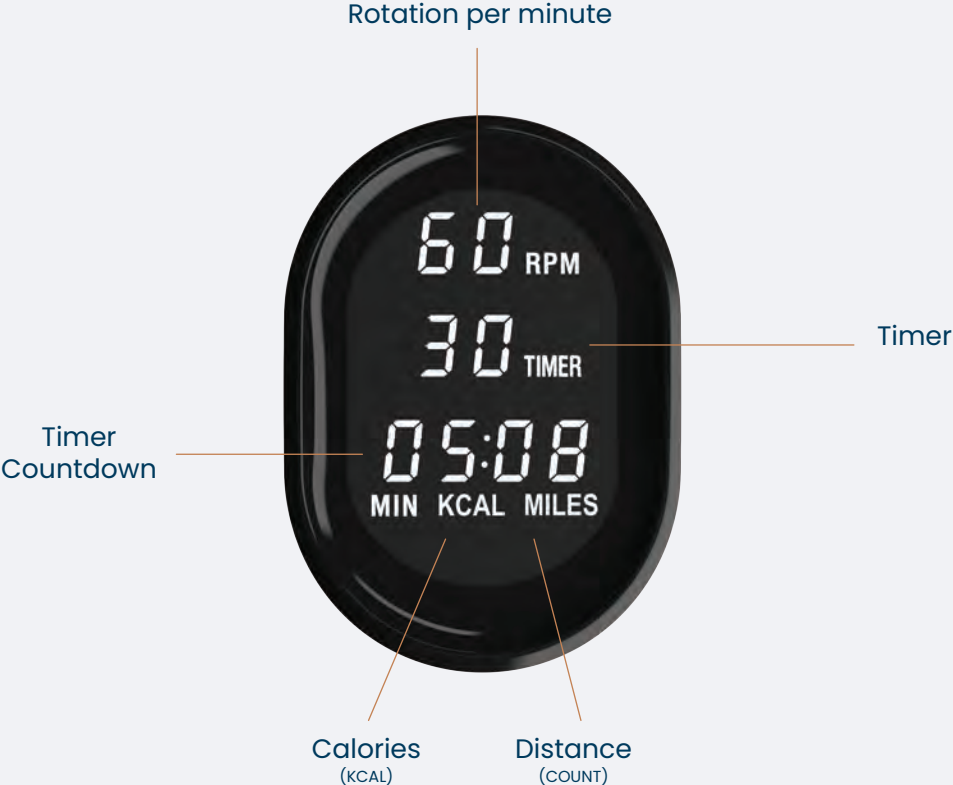


REMOTE

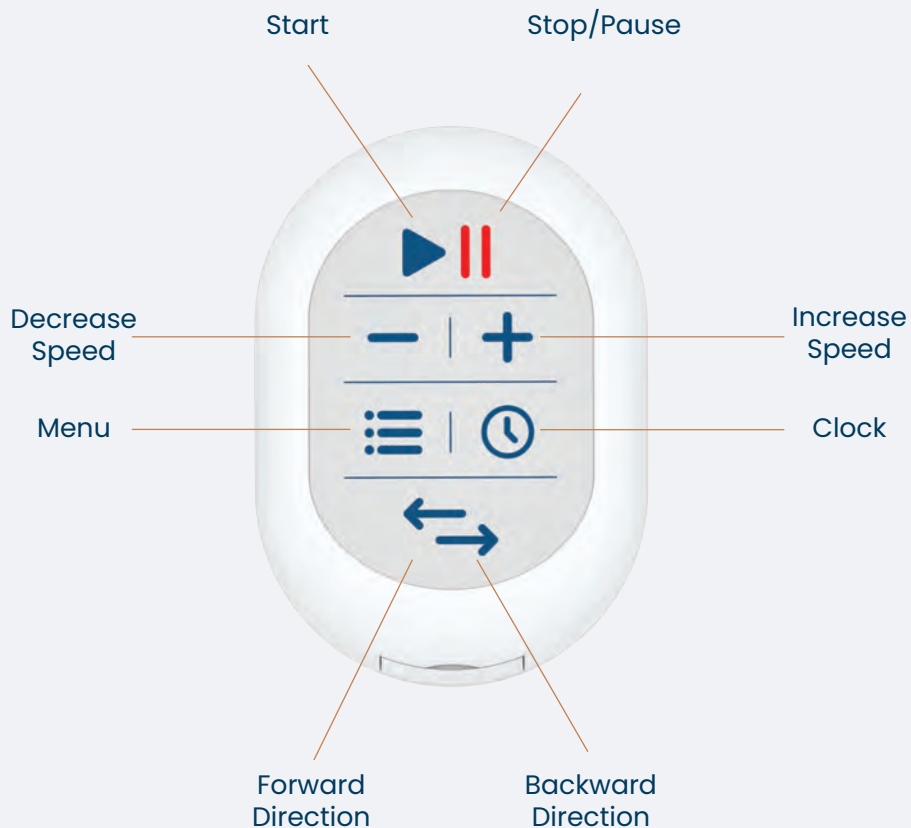


To use the PhysioPedal, screw the base onto the body by twisting the screws (x4).

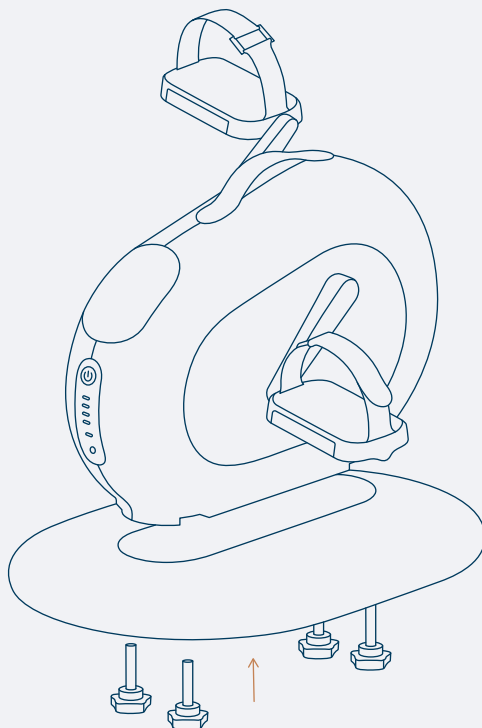

Attach the Pedal or Handles into the Cranks using The Quick Release Mechanism.

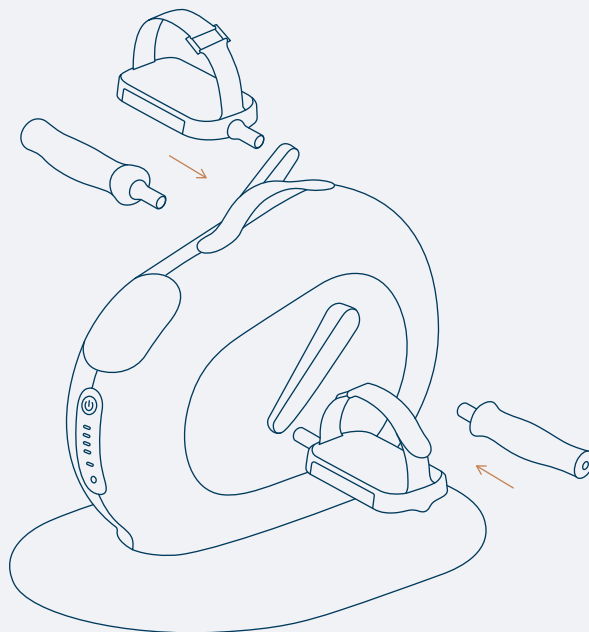


Insert the battery into the body by pushing it in until it clicks.

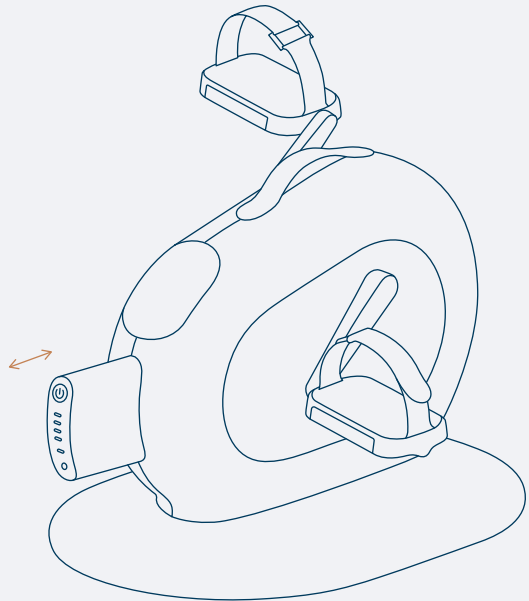


GETTING STARTED

To Power on the PhysioPedal, hold the power button down for 5 seconds until the battery lights up.

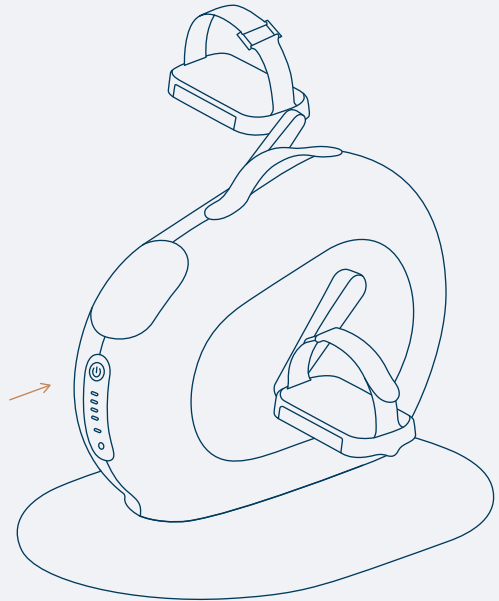


HOW TO USE

To exercise the upper body put the PhysioPedal on an elevated platform or table and put your hands on the grips.

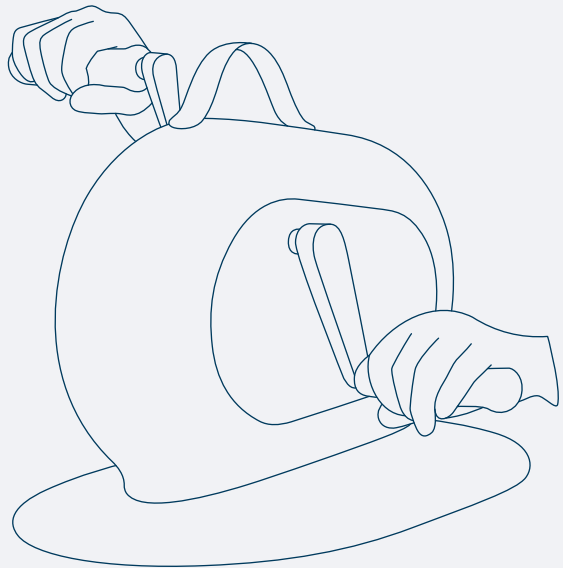

HOW TO USE

To use the PhysioPedal put both feet on the pedals.

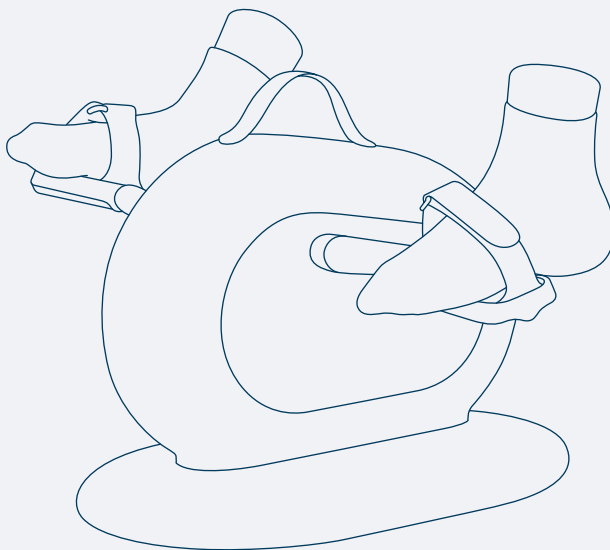

Press the START Icon to start the Pedals or Handles to move.




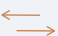
The unit will start a slow speed of 30 RPM, to increase the speed press the PLUS button. To decrease the speed press the MINUS button.



To set the timer, press the CLOCK Icon. The PhysioPedal will start at 10 minutes. Each press will increase the time by 10 minutes up to 90 minutes total.



To stop the Unit press the STOP/PAUSE Icon.

If the remote control stops working, you may need to reconnect the remote control to your PhysioPedal. Simultaneously press and hold the Power button  on the PhysioPedal and the Arrow button  on the remote control.

Make sure the power is OFF and the battery charger is not attached. Use a damp, clean cloth and gently wipe your PhysioPedal. Store in a clean, cool, dry place, not in direct sunlight, when not in use. Any other servicing should be performed by an authorized service representative.

If within three years from the date of purchase, this device fails due to a defect in material or workmanship, Nobol Inc. will repair or replace the device or necessary components, free of charge.

This warranty excludes:

- (a) Damage caused by accident, abuse, mishandling, or transport
- (b) Units taken out of the US
- (c) Units subjected to unauthorized repair
- (d) Damage exceeding the cost of the product
- (e) Deterioration of the delivered product resulting from abnormal storage and/or safeguarding conditions on the client's premises
- (f) Products taken apart or repaired by anyone other than an authorized service representative
- (g) Failure to provide the dated proof of purchase and serial number on your device
- (h) Products purchased from an unauthorized retailer
- (i) Please do not try and repair the unit yourself, as this will void your warranty.

Some states and countries do not allow a limitation of damages, so the foregoing limitation may not apply to you. This warranty guarantees specific legal rights, but other rights may vary between countries and from state to state. This warranty service is available through the operating policies and procedures of Nobol Inc., and all Nobol international distributors. If you are a client, please contact the appropriate Nobol business entity dealer, or retailer, from whom you directly purchased the device for warranty and return authorization procedures. This warranty is initiated and executed by Nobol Inc.

by **Nobel**[®]

Nobel[®]
Choose Health