




<p><b>Product</b> <b>Item number</b></p>	<p>MCT 3 change trolley, electrically height adjustable R9662721000</p>
<p><b>Color</b>  <b>Product description</b></p>	<p>EAN:5708590390000</p>  <p>White</p> <p>Pressalit offers a variety of mobile and wall-hung, height-adjustable and fixed height shower and changing tables of high quality. They are designed to provide a stable platform for showering and changing adults and children with reduced mobility. These tables can be used in schools, institutions, hospitals and private homes and are all designed to be durable, comfortable and easy to clean. The design and flexibility of the tables ensures optimal working conditions for the carer as well as a high level of comfort and safety for the user.</p> <p>MCT 3 change trolley, 27.6" x 82.7", two fold-down sides. Electrically height adjustable with wired remote control. 19.7" lifting range, from 17.7" to 37.4".</p> <p>Mattress to be ordered separately.</p>
<p><b>Adjustability</b> <b>Adjustment height</b>  <b>Product capacity</b>  <b>Speed up</b>  <b>Noise level</b></p>	<p>Height adjustable</p> <p>19,7 inch</p> <p>Maximum user weight (MUM): 441 lbs. Maximum load (SWL): 452 lbs. <i>The maximum load is to be understood as the total safe working load of the product, including both the weight of the user, and the weight of the mattress</i></p> <p><i>This product has been tested statically with a load of 1130 lbs as the IEC 60601-1 standard requires testing with a safety factor 2.5.</i></p> <p>Speed up: 0,67 inch/s Speed down: 0,67 inch/s</p> <p>50 dB</p>
<p><b>Materials</b></p>	<p>Plastics, compact laminate and powder-lacquered steel.</p>



<b>Product</b>	MCT 3 change trolley, electrically height adjustable
<b>Item number</b>	R9662721000
<b>Ergonomy</b>	<ul style="list-style-type: none"><li>• The ergonomically designed safety rails have a practical lower central section, which provides better access and work environment for the carer while transferring, changing and dressing the user.</li><li>• One hand operated safety rails allows the carer to maintain eye contact as well as physical contact with the user during the whole process, the foldable safety rails also provide easy access from all four sides of the table.</li><li>• The cut-out grips in the safety rails enable users to assist when being turned, which motivates physical training.</li><li>• Easy to clean, smooth contact surfaces which do not collect dirt and bacteria.</li></ul>
<b>Features</b>	<ul style="list-style-type: none"><li>• The trolley can be moved from room to room for easy transport and multiple users.</li><li>• Intuitive, easy-to-use drop down rails with one hand operation.</li></ul>
<b>Mains connection (V)</b>	120
<b>Actuator (motor)</b>	LA40 (Linak) IPX6
<b>Control box</b>	CA30 (Linak) IPX6
<b>Wired remote control</b>	HB30 (Linak) IPX6
<b>Power consumption (maximum load)</b>	Power consumption by maximum load: 330 W Power consumption by standby: 0,5 W
<b>Warranty</b>	Pressalit makes a commitment to remedy material and manufacturing defects for a period of three years from the invoicing date. NOTICE If the product is placed in an atmosphere heavy in salt or chloride, the warranty will be reduced from the normal 3 years to 1 year. Salt and chloride have a corrosive effect on the products components and will cause a decrease in the life expectancy of the product.
<b>Cleaning</b>	<p>This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p>
<b>Declaration of conformity</b>	This product is solely intended for the US and Canada, and has no EU declaration of conformity.

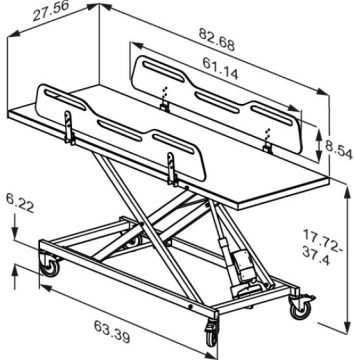

PRESSALIT®

<b>Product</b>	MCT 3 change trolley, electrically height adjustable
<b>Item number</b>	R9662721000