

# Performs - in Tune with You









# F&P ESON NASAL MASK


F&P Eson is designed to perform in tune with you and your needs.

Three small, streamlined mask components work in harmony for greater performance, stability and auto-adjusting comfort. This compact nasal mask is packed with technological innovations orchestrated to make SleepLife easier.





F&P Eson is very 'easy on you' from bag to bed. It has been designed to perform in tune with every aspect of SleepLife:\* F&P Eson maximizes no-fuss comfort and simplicity by minimizing manual adjustments, parts and noise.

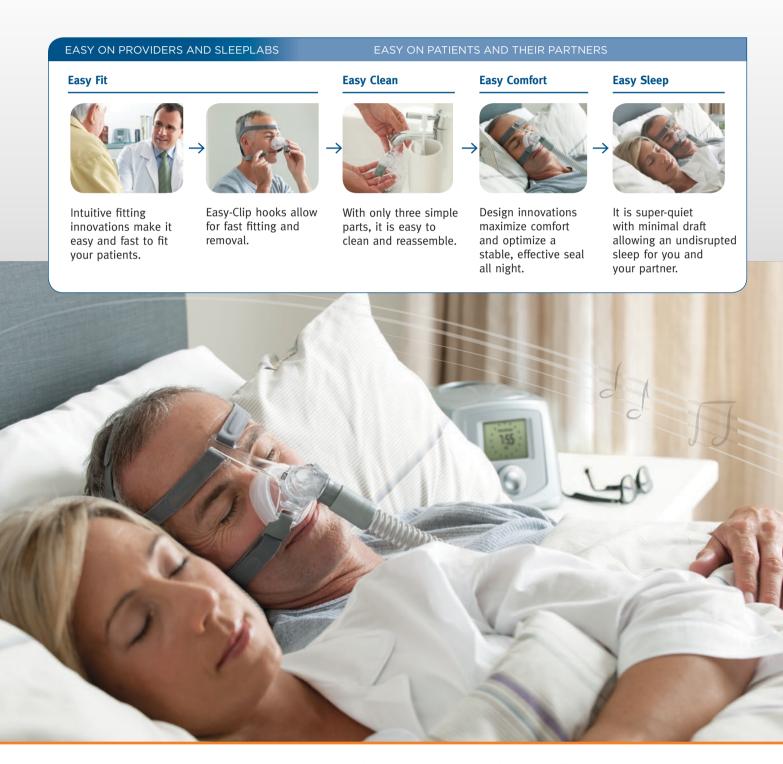

F&P Eson's three simple components, the **RollFit**<sup>™</sup> **Seal**, **ErgoFit Headgear** and **Easy Frame**, work in harmony to deliver the comfort, seal and easy use that Fisher & Paykel Healthcare masks are known for.



### 1. RollFit Seal

*Easy auto-adjusting fit that optimizes an effective, comfortable seal* 

As the name suggests, the onepiece seal 'rolls' back and forth on the bridge of the nose to adjust automatically. In doing so, the **RollFit Technology** minimizes pressure on the bridge of the nose.

The **RollFit Seal** incorporates **Integrated Stability Panels** for optimal sealing performance and stability.

### 2. ErgoFit Headgear

3a

Easy to use, fit, adjust and clean

2

This breathable, ergonomically designed headgear has been designed to self-locate high on the rear of the head allowing for maximum head movement (sideways, up and down) without mask dislodgement.

The **ErgoFit Headgear** also incorporates **Embedded VELCRO®\* Tabs** (2a) with **Easy-Clip hooks** (2b) for fast and easy fitting and removal. The adjustable **Crown Strap** (2c) enhances mask stability.

\* VELCRO® is a registered trademark of Velcro Industries B.V.

### 3. Easy Frame

2b

Easy to wear, see and sleep

2a

This low-profile frame is stable, durable and small, and ensures a clear line of sight.

The one frame fits all three seal sizes and has an **'Easy-Clip'** frame attachment for effortless assembly after cleaning.

The **Easy Frame** comes with the **Q Cover** (3a) for super-quiet use and minimal partner draft.

The Ball and Socket Elbow rotates for freedom of movement while minimizing drag from the CPAP tube.

F&P Eson performs in tune with what patients and their healthcare providers want: **Performance** that makes adherence **Easy**.

# For more information, please contact your local Fisher & Paykel Healthcare representative

### Manufacturer

Fisher & Paykel Healthcare Ltd 15 Maurice Paykel Place East Tamaki, Auckland 2013

PO Box 14 348 Panmure Auckland 1741 New Zealand

+64 9 574 0100 Tel· Fax: +64 9 574 0158 Email: info@fphcare.co.nz Web: www.fphcare.com

#### Australia

Fisher & Paykel Healthcare Pty Limited 36-40 New Street, PO Box 167 Ringwood, Melbourne Victoria 3134, Australia

Tel: +61 3 9879 5022 Fax: +61 3 9879 5232

#### Austria Tel: 0800 29 31 23

Fax: 0800 29 31 22

## Benelux

Tel: +31 40 216 3555 Fax: +31 40 216 3554

## China

Tel: +86 20 3205 3486 Fax: +86 20 3205 2132

# France

Tel: +33 1 6446 5201 Fax: +33 1 6446 5221

### Germany

Tel: +49 7181 98599 0 Fax: +49 7181 98599 66 India Tel: +91 80 4284 4000 Fax: +91 80 4123 6044

**Irish Republic** Tel: 1800 409 011

Italy Tel: +39 06 7839 2939 Fax: +39 06 7814 7709

Japan Tel: +81 3 3661 7205 Fax: +81 3 3661 7206

Northern Ireland Tel: 0800 132 189

Spain Tel: +34 902 013 346 Fax: +34 902 013 379

#### Sweden Tel: +46 8 564 76 680 Fax: +46 8 36 63 10

Switzerland Tel: 0800 83 47 63 Fax: 0800 83 47 54

### Taiwan

Tel: +886 2 8751 1739 Fax: +886 2 8751 5625

#### Turkey

Fisher Paykel Sağlık Ürünleri Ticaret Limited Şirketi, P.O. Box 06371 Ostim, Alinteri Bulvari 1161/1 Sokak No. 12-14, Ankara, Turkey

Tel: +90 312 354 34 12 Fax: +90 312 354 31 01

### UK EC REP

Fisher & Paykel Healthcare Ltd Unit 16, Cordwallis Park Clivemont Road, Maidenhead Berkshire SL6 7BU, UK

Tel: +44 1628 626 136 Fax: +44 1628 626 146

#### USA/Canada

Tel: +1 800 446 3908 +1 949 453 4000 or Fax: +1 949 453 4001

フィッシャー&パイケルヘルスケア株式会社 〒103-0025 東京都中央区日本橋茅場町3-8-10 リベラ茅場町5F FAX. 03-3661-7206 ご注文書専用 FAX. 0800-100-7207 TEL. 03-3661-7205 (代表)

| Mask Reorder Codes |                      | Spare Parts Reorder Codes |                                    |          |                                    |
|--------------------|----------------------|---------------------------|------------------------------------|----------|------------------------------------|
| 400449             | F&P Eson Mask Small  | 400HC567                  | F&P Eson Headgear - Small          | 400HC228 | F&P Eson Diffuser (10pk) and Cover |
| 400450             | F&P Eson Mask Medium | 400HC568                  | F&P Eson Headgear - Medium/Large   | 400HC570 | F&P Eson Swivel                    |
| 400451             | F&P Eson Mask Large  | 400HC569                  | F&P Eson Headgear - Clips & Buckle | 400HC571 | F&P Eson Mask No H/G - Small       |
|                    |                      | 400HC126                  | F&P Eson Seal - Small              | 400HC572 | F&P Eson Mask No H/G - Medium      |
|                    |                      | 400HC127                  | F&P Eson Seal - Medium             | 400HC573 | F&P Eson Mask No H/G - Large       |
|                    |                      | 400HC128                  | F&P Eson Seal - Large              |          |                                    |

