

### **HEALTHCARE SEATING**

**Patient Comfort** Caregiver Functionality





# **LUMEX**<sup>®</sup>

Since 1946, LUMEX® has manufactured quality healthcare products for use in hospitals, clinics, extended care facilities and the home. Lumex has been the innovator in healthcare seating solutions for use in a wide variety of healthcare applications:

- Dialysis and Nocturnal Dialysis
- Oncology and IV Therapies
- Same Day Surgery
- Post Operative Care
- Post Anesthesia Recovery (PACU)
- Patient Room Seating
- Bariatric Care
- Blood Collection
- Extended Care
- Sub-Acute Care

LUMEX Healthcare Seating can help liberate almost anyone from the confines of a bed and help make them more comfortable for longer periods of time. Lumex Clinical Care Recliners makes the care, treatment and recovery more comfortable for patients and caregivers alike. For extended care residents our recliners can be modified with options to meet special needs. Lumex Clinical Care Recliners are more versatile and mobile than a bed and can be used throughout a healthcare facility.

## All Lumex Healthcare Seating is protected by our Limited Wraparound Warranty\*:

Welded Frame: Lifetime<sup>1</sup> Warranty

Mechanical Components: Three (3) Years

Upholstered Components: One (1) Year

Casters: One (1) Year

Heat and Massage Option: Two (2) Years

Non-Welded Frames: Three (3) Years

\*Subject to Warranty Statement Limitations

¹(Defined as Ten (10) Years)

For the most update information on all GF Health Products, Inc., go to *www.grahamfield.com*.



The Original Lumex "Wheeled Chair Table" circa 1986



A "Beauty Shot" from a 1986 Brochure



Cutaway Drawing of a Lumex Six Position Recliner circa 1982



The Original Lumex 573 Three Position and Lumex 576 Six Position Recliners



The Original Lumex Six Position Recliner circa1984

### **LUMEX CLINICAL CARE RECLINERS**

Lumex Clinical Care Recliners offer healthcare facilities specialized seating that benefits both patients and caregivers. Designed and manufactured in the USA, Lumex Clinical Care Recliners feature easily replaceable components that increases the lifespan of the recliner and reduces their cost in use over time.

### **Features:** High Density Foam

For maximum patient comfort during any time intensive healthcare procedure all Lumex Clinical Care Recliners feature high-density foam in the seat, back and armrests. All upholstery components are easily replaced should they become worn or damaged.

**Healthcare Grade Vinyl Upholstery** 

For beauty, long life and ease of cleaning and disinfection, Lumex Clinical Care Recliners are available in a beautiful palette of colors / fabrics.
The Lumex Clinical Care Recliner Fabric selection can address the concerns of end users, facility management and interior designers.

### **Standard Vinyl Collection:**

- Patriot Plus: The Vinyl Upholstery "Standard" designed specifically for healthcare applications in a beautiful palette of colors
- Renaissance: Durable vinyl with a subtle pattern
- Savoy: Same features as Patriot Plus Vinyl in a softer "hand"




SilverGuard Vinyl Collection: The SilverGuard Collection features SILVERGUARD® Silver Ion Technology. SILVERGUARD is a natural, environmentally friendly bacteria-fighting defense to inhibit odor and stain causing bacteria. SILVERGUARD technology is EPA registered.

- Dolce: Soft, leather-like urethane based material for maximum patient comfort
- Guardian: Minimal texture, exclusive top coat and Silver Ion Technology make this our toughest upholstery offering
- Silvertex: A beautifully reflective fabric featuring Silver Ion Technology.

### **How SILVERGUARD® works:**

Silver lons destroy micro-organisms by attaching to the bacteria cell walls to disrupt protein metabolism within the cells. These are slowly released from their inorganic matrix in the

presence of moisture. Simply stated-wipe down the upholstery with a damp cloth to activate the Silver lons.

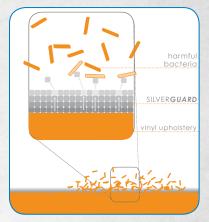

### TENTE® Casters

The mobility of any recliner is paramount in a healthcare setting. Moving patients



from a procedure room to a recovery area or repositioning a recliner in a tight dialysis station demands a high quality caster. Every Lumex Clinical Care Recliner utilizes TENTE brand casters, world renowned for their quality and durability. TENTE casters are quiet, non-marking and virtually maintenance free.

Lumex Clinical Care Recliner features three swivel, total lock casters identified by a black locking tab. The fourth caster, identified by a green locking tab, is a directional locking caster that when actuated, allows the recliner to track straight when being pushed down long hallways.

### Fold Down Side Tables

To allow for easier patient preparation, Lumex Clinical Care Recliners feature two plastic fold-down side tables. Each side table features a smooth surface with no sharp edges to



catch dirt making cleaning and disinfection easier. Side tables sit evenly with armrest in the up position and when not in use, the tables fold flush within the sides of the recliner to more easily fit through doorways.

### **Heat and Massage Function**

During prolonged treatment, patients tend to feel cold, particularly while undergoing dialysis. Lumex offers an integrated heat and massage unit on select clinical care recliners. The heat and massage function is easily accessible by the patient or caregiver and can be operated independently of each other to give patients a choice in managing either heat, massage or both. Heat and massage is available on the following models:

- FR577RG
- FR566G / FR566DG
- FR587W



An easy-to-use switch independently controls both the Heat and Massage features. Two heating elements, one in the seat and one in the back, provide gentle heat for the patient. Three heat settings allow patients to control just how much warmth they need. The heat feature has an automatic shut-off to help ensure that the patient is always in control of their comfort.

### **CA TB-133 Flammability Testing**

All Lumex Clinical Care Recliners meet the stringent California Technical Bulletin 133 fire rating. Whether your state requires upholstered furniture to meet this stringent test or not, knowing all Lumex Clinical Care Recliners meet these test parameters will give you peace of mind. CA TB-133 compatibility is available as an option on most Lumex Long Term Care Recliners.

Features a rear mounted Trendelenburg release lever and specially designed back mounted arm pads to support the patient's arms while they are reclined. This recliner is a cost-effective choice for all clinical care applications. Wider and softer armrests for patient comfort feature wall-saving construction to protect walls and helps to prevent damage to the armrest upholstery.

Recline Positions: Features three positions for patient comfort; the fourth position is the Trendelenburg (or shock) position that must be activated by a caregiver. The Trendelenburg Position cannot be activated by a patient while they are sitting in the recliner.

- First or Upright Position: Ideal for patient transport and to initiate patient treatment / therapy.
- Second or First Recline Position: Ideal for patient treatment, rest, TV and reading.
- Third or Full Recline Position: Ideal for patient treatment, rest and reading.
- Trendelenburg Position (Caregiver activated only)

### Unique Features of the Series FR577RG and FR577RGH:

- Higher Floor-to-Back Height: Enables the caregiver to raise the patient up with less effort, reducing the chances for a back injury.
- Trendelenburg Position: Easily actuated by the Caregiver from the back of the recliner via a trigger release mechanism mounted on the back push bar.
- Side Panels: Easily removed side panels enable access to the recliner for cleaning and maintenance.
- Heat and Massage Option: Available
- Testing: California Technical Bulletin CA TB-133 Flammability Test
- Weight Capacity: 290lb Evenly Distributed





A specially contoured seat and back for the ultimate in patient comfort. The extra wide back design comfortably supports the patient from upright through each recline position. Foot activated Trendelenburg pedals located on either side of the recliner allows the caregiver to stay in eye contact with the patient while placing the patient into the Trendelenburg Position.

Recline Positions: Features three positions for patient comfort; the fourth position is the Trendelenburg (or shock) position that must be activated by a caregiver. The Trendelenburg Position cannot be activated by a patient while they are sitting in the recliner.

- First or Upright Position: Ideal for patient transport and to initiate patient treatment / therapy.
- Second or First Recline Position: Ideal for patient treatment, rest, TV and reading.
- 3. **Third or Full Recline Position:** Ideal for patient treatment, rest and reading.
- 4. Trendelenburg Position (Caregiver activated only)




- Trendelenburg Position: Easily actuated from either side of recliner by conveniently located foot pedals. Allows the caregiver to stay in constant eye contact with the patient while repositioning the patient.
- Side Panels: Easily removed side panels enable access to the recliner for cleaning and maintenance.
- Heat and Massage Option: Available
- Testing: California Technical Bulletin CA TB-133
   Flammability Test
- Weight Capacity: 350lb Evenly Distributed



With the same features as the Series FR566G, the Series FR566DG is made specifically for patients who cannot easily Stand, Pivot and Sit. This recliner features an easy-to-release drop-arm mechanism on both arms that positions the arm level with the seat to allow for easier lateral transfer from a bed, wheelchair or stretcher.

### Unique Features of the Series FR566DG and FR566DGH:

- Drop Arm Feature: Allows both arms to drop down evenly with the seat allowing for easier and safer lateral transfers for both the patient and caregiver.
- Trendelenburg Position: Easily actuated from either side of recliner by conveniently located foot pedals. Allows the caregiver to stay in constant eye contact with the patient while repositioning the patient.
- Side Panels: Easily removed side panels enable access to the recliner for cleaning and maintenance.
- Heat and Massage Option: Available

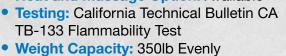





Featuring an increased weight capacity of 450 pounds (evenly distributed) and a wider seat, the Lumex Deluxe Clinical Care Recliner-Wide is ideal for bariatric patients and clinics. The extra-wide 24" seat width and greater weight bearing capacity enables the recliner to comfortably sit and position patients up to 450 pounds (evenly distributed). The recliner mechanism features dual gas cylinders that give smooth power assistance for back position adjustment and Trendelenburg activation.

Recline Positions: Features three positions for patient comfort; the fourth position is the Trendelenburg (or shock) position that must be activated by a caregiver. The Trendelenburg Position cannot be activated by a patient while they are sitting in the recliner.

- First or Upright Position: Ideal for patient transport and to initiate patient treatment / therapy.
- Second or First Recline Position: Ideal for patient treatment, rest, TV and reading.
- Third or Full Recline Position: Ideal for patient treatment, rest and reading.
- Trendelenburg Position: (Caregiver activated only)

### Unique Features of the Series FR587W and FR587WH:

- Trendelenburg Position: Easily actuated by the Caregiver from the back of the recliner via a trigger mechanism mounted on the back push bar.
- Side Panels: Easily removed side panels enable access to the recliner for cleaning and maintenance.
- Heat and Massage Option: Available
- Testing: California Technical Bulletin CA TB-133 Flammability Test
- Weight Capacity: 450 lb Evenly Distributed




Combining beauty and full clinical functionality, this recliner has been specifically designed for use in a wide variety of healthcare environments. Featuring an independent legrest function and infinite positioning of the back, patients will appreciate both the comfort and functionality of the recliner. The Lumex Ortho-Biotic II Recliner is proudly Made in the USA.

Recline Positions: Features independent legrest operation and infinitely adjustable back within a range. Trendelenburg Position is easily actuated by side mounted foot pedal allowing the caregiver to stay in constant contact with the patient while repositioning the patient.

- 1. First or Upright Position: Ideal for patient transport and to initiate patient treatment / therapy.
- 2. First or Upright Position with Legrest Extended: Ideal for patient transport and rest.
- 3. Infinitely Adjustable Back Position
- 4. Trendelenburg Position: (Caregiver activated only)

### **Unique Features of the Series FR597G and FR597P:**

- Legrest-Recline Lever: Conveniently located on the outside of the recliner, the lever enables the patient or caregiver to easily operate the legrest and adjust the back angle of the recliner.
- Wood Capped or Fully Upholstered Arms:
   Beautifully finished wooden armrests add a touch of elegance and make the Lumex Ortho-Biotic II perfect for post-operative care and as a patient room recliner. Also available with fully upholstered arms, Series FR597P.
- Self-Storing Footrest: For patient comfort and enhanced safety during patient transport, a pullout footrest is standard. The footrest is easily extended when needed and can be safely stored out of the way under the recliner when not in use.
- Universal Headrest: The Universal Headrest can easily be adjusted up and down on the back of the recliner. Unique shape allows the pillow to be used behind the patient's head, neck or back.
- Composite Foam Seat and Back: The seat and back are constructed of a high density foam base with a 1" topper of visco-elastic foam.

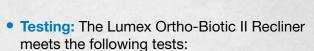

- California Technical Bulletin CA TB-133 Flammability Test
- ANSI/BIFMA X5.4-2013 Sections 11,12,14, & 18
- Weight Capacity: 350lb Evenly Distributed
- Factory Installed Accessories:
- Side Table: FR597GTBLMNT
- IV Pole Mounting Bracket: FR597GIVMNT
- IV Pole: 2101A



### **Lumex Clinical Care Recliner Accessories**

Universal Headrest Pillow: Item FR5659US

The Universal
Headrest Pillow's
unique shape
features an angled
design that offers the
patient a different level
of support depending
on whether the small contour or
large contour is placed under the
patient's head.

- Easily positioned for individualized patient comfort
- Easy-on / easy-off hook and loop attachment
- Can be positioned anywhere along back of recliner
- Fits All Lumex Recliners

### **Headrest Cover**

Headrest covers are the ideal way to protect and extend the life of your Lumex Recliners. They help protect the vinyl from hair oil, hair spray, coloring agents, etc., the leading causes of premature aging of the vinyl. Available in matching colors, the headrest covers feature an easy-on / easy-off design for cleaning or replacement. Available for the following series:

- FR566G / FR566GH
- FR566DG / FR566DGH
- FR577RG / FR577RGH
- FR587W / FR587WH

### **How to Order:**

• HRC566: Fits Series FR566G &

FR566DG Recliners

HRC577: Fits Series FR577RG Recliners

HRC587W: Fits Series FR587W Recliners



### IV Pole Attachment: Item #: 5782RG

Offers a tip-free means of administering enteral feedings or IV therapy. For use with Lumex Bed Socket IV Pole (Model 2101A); sold separately

### **Fits Series:**

FR565WG / FR566G FR577RG / FR587W Does NOT Fit Series FR566DG Recliners



### Bed Socket IV Pole: Item: 2101A (1/CTN) and 2100A (3/CTN)

Locking collet quickly and securely adjusts height from 35" to 65". Internal rubber bumpers ensure silent operation and reduce shock to hand and I.V. container when lowering telescoping tube. Constructed of durable, anodized aluminum tubing.

#### Fits:

- 5782RG
- FR597GIVMNT

### Leg Extensions: Item:

Fits all round tube frame designs.

- Adds 5" to height of recliner.
- Raises recliner by five Inches (5")

#### **Fits Series:**

- FR565WG
- FR566G
- FR577RG
- FR587W

NOTE: Does

**NOT Fit Series** 

FR566DG

Recliners



This Recliner can be a great assist to caregivers caring for Bariatric Patients. Featuring a heavy duty welded steel frame and heavy duty actuators, this recliner can more safely position patients of up to 700 pounds (when evenly distributed).

Recline Positions: Features infinite positions with a range for patient comfort via a four-function hand pendent. Legrest and backrest can be raised and lowered independently, giving the caregiver complete control over positioning. Powered unit is equipped with a rechargeable battery pack for emergency operation during a power failure.

### **Unique Features of the Series FR588W:**

 Trendelenburg Position: Quick-release trigger mechanism on push handle allows chair to swiftly move into Trendelenburg Position.

 Drop Arm Feature: Drop arm mechanism allows both arms to drop down allowing for easier and safer lateral transfers for both the patient and caregiver.

 Patient Footplate: Enables easier entry and exit into the recliner and safely positions the patient's feet during transport.

 Wall Saver Armrest Caps: Increases the durability of the armrests and protects walls from marring.

 Fold Down Side Tables: Two, durable, high-impact side tables, easy to actuate by caregiver, sit flush with recliner arm to conveniently hold prep materials; side table surface is smooth and features tapered edges for easy cleaning.

- Heavy Duty Casters: Four swivel casters facilitate transport and maneuverability in confined areas.
- Testing: California Technical Bulletin CA TB-133 Flammability Test
- Weight Capacity: 700 lb Evenly Distributed

### **Available Colors:**

Blue Ridge: FR588W427
 Doe Skin: FR588W851
 Jade: FR588W857
 Rosewood: FR588W863



Ideal for bariatric care or patient requiring an extra-wide 24" wide seat, this recliner features a weight capacity of 400 lb evenly distributed. Features infinite positioning within a range and a unique position lock that can stabilize the position of the recliner for patient comfort, and safety. The patented 'link and pin' mechanism design allows for smooth position change. Ideal for use as a Patient Room Recliner, Post Operative Care Recliner, or for Oncology and Blood Collection.

Recline Positions: The Deluxe Preferred Care Recliner can attain almost infinite positioning within a range for patient comfort.

### **Unique Features of the Series FR565WG:**

- Wall Saver Armrest Caps: Increases the durability of the armrests and protects walls from marring.
- Quick-Release Upholstery: Upholstery can be quickly and easily removed for easy cleaning and replacement.
- Ultra-Flex Suspension Bands: By simply removing the back upholstery from the recliner back, the Ultra-Flex Suspension Bands can be adjusted to meet individual patient comfort requirements.
- Testing: California Technical Bulletin CA TB133 Flammability Test
- Weight Capacity: 400 lb Evenly Distributed

### **Accessor:**

Front Tray Table: Model 5644WG



A truly multi-faceted product, the Lumex Preferred Care® Bed-Recliner has been designed for a wide variety of healthcare applications:

- Dialysis
- Nocturnal Dialysis
- Sub-Acute Care
- Oncology / IV Therapies
- Rehabilitation
- Emergency Room

The Preferred Care
Bed-Recliner is
ideal for use
during any
intensive
health-care
procedure or
when a versatile
combination
of a bed and



recliner is called for. Unlike recliners that are converted into a semi-flat sleeping surface, the PCB-R offers the comfort of a true bed with an extremely comfortable sleep surface. From a bed, the PCB-R can attain all the positions of a recliner with an almost infinite number of seated and reclined positions.

### **Unique Features of the AR981642:**

- Low position of 7.0" to a Deck Height to 30" Deck Height (without Sleep Surface and Overlay). This wide range of motion allows for easier care by the caregiver in the full up position and helps prevent patient injuries from falls while in the low position.
- Cardiac Position: Is primarily designed for residents with heart and respiratory problems. Helps to elevate the patient from the waist up, which provides relief to the lungs, increases circulation, and assists the patient in coming from a fully prone position to a sitting position without undue strain.



**Trendelenburg Position:** Full or Reverse Trendelenburg Positions and Auto CPR Position

- All functions can be "locked-out" and managed by the caregiver for patients with cognitive challenges.
- Sleep Surface and Overlay: A unique two-piece design ensures the greatest support and comfort for the patient in any prone or reclined position.

Sleep Surface: Constructed with 3.5" High Density Base Foam and a 1.5" Top Wrap Comfort Layer, the sleep surface is affixed to the moving sections of the bed to ensure support in every position. Upholstered with Silvertex Silver Ion Vinyl.

Overlay: Constructed with a 1.5" high density base foam and a 1.5" Visco-Elastic top foam for the ultimate in patient comfort. The Overlay is upholstered with a top cover of fluid-resistant five way stretch fabric and a fluid resistant non-skid bottom cover. Elastic attachment straps at each corner allow the Overlay to move in conjunction with Sleep Surface. Can be easily removed for cleaning.

**Testing:** The Lumex Preferred Care Bed-Recliner is built to conform to UL 60601-1 and IEC601-2-38; CAN/CSA C22.2 No. 601.1. Multiple Patents Pending

Weight Capacity: (including resident, support surface & accessories): 600 lb Evenly Distributed.

| Specifications                                   | Dimensions |
|--|------------|
| Mattress Deck Length                             | 80"        |
| Mattress Deck Width                              | 35"        |
| Overall Bed Length<br>Without Wall Saver         | 85.5"      |
| Overall Bed Width                                | 41.0"      |
| Maximum Deck Height                              | 30"        |
| Minimum Deck Height                              | 7"         |
| Maximum Deck Height w/Sleep Surface & Overlay    | 37.5"      |
| Minimum Deck Height<br>w/Sleep Surface & Overlay | 14.5"      |
|  |            |

Safe Workload (including resident, support surface & accessories): 600 lb Max Weight capacity
Evenly Distributed.

The Lumex® Preferred Care® Bed-Recliner is built to conform to UL 60601-1

The standard of the healthcare industry and the first recliner specifically designed to improve the comfort and mobility of residents in extended care facilities.

Offers numerous therapeutic benefits including position change that can aid in circulation and ease discomfort from respiratory, heart and other ailments. Position stabilization bar allows the caregiver only to secure the recliner in any one of the three positions.

For safer, easier transport, the Lumex Three Position Recliner features five-inch rear swivel casters with wheel locks and five-inch fixed front wheels. The molded activity table can be fixed in position or left unsecured for easy exit from the recliner. Stores at side of recliner when not in use. Durable welded steel frame.

#### **Recline Positions:**

- First or Upright Position: Ideal for patient transport and to initiate patient treatment / therapy.
- Second or First Recline Position: Ideal for patient treatment, rest, TV and reading.
- Third or Full Recline Position: Ideal for patient treatment, rest and reading.
- Available in standard upholstery or CA TB-133 rated upholstery.

Weight Capacity: 250lb Evenly Distributed

### **Caster Upgrade:**

 Rear caster upgrade – Tente casters feature total swivel lock. Easily installed. 5740TCSTR (2 required)

### Available Colors: Standard Upholstery

Vintage Gold: 574G401 Warm Taupe: 574G409 Blue Ridge: 574G427 Imperial Blue: 574G432 Royal Blue: 574G454 Doe Skin: 574G851 Jade: 574G857 Rosewood: 574G863

### **Available Colors:**

### **CA TB-133 Upholstery**

Vintage Gold: FR574G401 Warm Taupe: FR574G409 Blue Ridge: FR574G427 Imperial Blue: FR574G432 Royal Blue: FR574G454 Doe Skin: FR574G851 Jade: FR574G857 Rosewood: FR574G863





### Standard features found on all Lumex Preferred Care Recliners:

- Reduced armrest and seat height enable the recliner to fit under dining / activity tables and allow residents to "scoot" more easily.
- Two swivel casters with two rear locking casters enable the Preferred Care Recliner to be easily maneuvered in small spaces and track straight when pushed down hallways.
- A unique Ultra-Flex suspension system in the seat and back adjusts to each person, cradling the resident from head to heel, greatly reducing the risk of "bottoming out".
- By simply removing the back upholstery from the recliner back, the Ultra-Flex Suspension Bands can be adjusted to meet individual patient comfort requirements.
- Extra-wide padded armrests and side panels maximize comfort.
- Quick-release upholstery allows for easy cleaning, removal and replacement.
- A variety of accessories is available to customize the recliners for specific resident needs.
- CA TB-133 Rated Upholstery
- Weight Capacity: 250lb Evenly Distributed

### Lumex Preferred Care Recliner-Series FR565G

### **Available Colors:**

CA TB-133 Upholstery

Blue Ridge: FR565G427Jade: FR565G857Rosewood: FR565G863



All of the features and benefits of the Preferred Care Recliner plus a unique drop arm mechanism. The self-storing, height adjustable drop arm is simple for the caregiver to activate and enables a safer, easier lateral transfer for the resident and caregiver.

### **Available Colors:**

CA TB 133 Upholstery

Blue Ridge: FR565DG427Jade: FR565DG857Rosewood: FR565DG863

### Lumex Preferred Care Recliner with Tilt-in-Space Function: Series FR565TG

The Tilt-in-Space feature allows the caregiver to change the "tilt" of the seat and back for pressure relief and position change. Please note that this version of the recliner does not recline but utilizes a side mounted foot pedal mechanism to extend the legrest and "tilt" the position of the patient while maintaining a constant back to seat angle. In addition, the recliner features a self-storing, height adjustable drop arm that is simple for the caregiver to activate and enables a safer, easier lateral transfer for the caregiver and resident.

#### **Available Colors:**

CA TB 133 Upholstery

Blue Ridge: FR565TG427Rosewood: FR565TG863



### **Lumex Preferred Care Recliner Accessories**

Soft, comfortable positioning accessories. Helps maintain comfortable body alignment. Accessories can be easily installed without tools.

### Item FR56461 Head Support Standard Fits Series:

FR565G / FR565DG / FR565TG

#### **Available Colors:**

Blue Ridge: FR56461427
 Jade: FR56461857
 Rosewood: FR56461863

### Body Bolster: Item FR56471 NOTE: Body Bolster WILL NOT fit Series FR565WG Recliners Fits Series:

FR565G / FR565DG / FR565TG

#### **Available Colors:**

Blue Ridge: FR56471427
 Jade: FR56471857
 Rosewood: FR56471863

### Activity Tray Table: Item 5644G and 5644WG

Strong polypropylene tray table offers a supportive, stable work surface for activities. Includes receiving arms, storage bracket, and strap. Easily stored in the back of chair. Patent pending "lock / lock-out" key gives flexibility to the caregiver.

#### **Item 5644G Fits Series:**

FR565G / FR565DG / FR565TG

### **Item 5644WG Fits Series:**

FR565WG



### **Lumex Prefered Care Recliner Accessories**

### IV Pole Attachment: Item 5782RG

Offers a tip-free means of administering enteral feedings or IV therapy. For use with Lumex Bed Socket IV Pole (Item #: 2101A); sold separately.

#### **Fits Series:**

 5FR565G / FR565DG / FR565TG / FR565WG

### Bed Socket IV Pole: Item 2101A (1/CTN) and 2100A (3/CTN)

Locking collet quickly and securely adjusts height from 35" to 65". Internal rubber bumpers ensure silent operation and reduce shock to hand and I.V. container when lowering telescoping tube. Constructed of durable, anodized aluminum tubing.

#### Fits:

- 5782RG
- FR597GIVMNT



Enhances patient positioning and comfort and aids in the prevention of foot drop. Footrests tip up and out of way to allow access for easy patient transfer. Attaches without the use of tools.

### **Fits Series:**

 FR565G / FR565DG / FR565TG / FR565WG



### **Optional Side Bolsters: Item 599G**

Adds lateral support for residents who have weak upper body control. Easily mounts to back of recliner.

#### **Fits Series:**

574G / FR574G



### Foot-Drop Assembly: Item 5657G

Helps keep resident's feet positioned correctly to help avoid foot-drop. NOTE: Foot-Drop Assembly will not fit FR565WG.

### **Available Colors:**

Blue Ridge: 5657G427Jade: 5657G857Rosewood: 5657G863



The Everyday Hip Chair is ideal for people who have arthritis or for patients who have undergone hip or knee surgery. The 27" high seat helps reduce flexion, joint stress and fatigue while sitting and when entering / exiting the chair.

### **Features:**

- High-density seat and back foam makes sitting more comfortable for longer periods of time.
- Adjustable height footrest provides greater comfort while sitting.
- Footrest height is easily and securely adjusted by two spring-loaded adjustment knobs located on either side of frame.
- Padded full-length arms make sitting and standing easier.
- The anti-slip strip on the adjustable footrest provides added safety and is an assist during repositioning.
- Sturdy metal frame and padded seat and back are easy to clean and maintain.

### **Available Colors:**

Blue Ridge: GF4405427
 Jade: GF4405857
 Rosewood: GF4405863





**Adjustable Height for Patient Comfort** 



### **Stretchers and Procedure Chairs**

Hausted has been designing and manufacturing a unique line of patient handling systems with its Stretchers and Procedure Chairs to be used in various procedures from eye surgery to bariatric care, radiology and more.





### **Therapeutic Support Surfaces**

One of the most well known names in sleep surfaces, Simmons Clinical Care, offers mattresses, alternating pressure / low air loss systems, overlays, cushions and replacement stretcher pads. Simmons provides both static and dynamic systems to address pressure management from prevention to treatment.





As the inventor of the folding wheelchair in 1933, E&J<sup>™</sup> Wheelchairs continues to offer cost effective patient transport solutions for a wide variety of patient mobility needs.



|  |        |       |         |       |         |       | _        | .umex    | Healt   | hcare  | Seati    | ng Din | Lumex Healthcare Seating Dimensions | JIS    |          |         |            |          |         |       |         |       |        |
|--|--------|-------|---------|-------|---------|-------|----------|----------|---------|--------|----------|--------|-------------------------------------|--------|----------|---------|------------|----------|---------|-------|---------|-------|--------|
| DIMENSION  | FR566G | 96G   | FR566DG | sDG   | FR577RG | 'RG   | FR587W   | <b>A</b> | FR588W  | >      | FR597    |        | 574G / FR574G                       | 74G    | FR565G   |         | FR565TG    |          | FR565DG | ш     | FR565WG |       | GF4405 |
|  | inch.  | cm.   | inch.   | cm.   | inch.   | cm.   | inch.    | cm.      | inch.   | cm.    | inch.    | cm.    | inch.                               | cm. in | inch.    | cm. in  | inch. cm.  | n. inch. | . cm.   | inch. | h. cm.  | inch. | cm.    |
| Overall Height:                                  | 49.5   | 125.7 | 50.5    | 128.3 | 48      | 121.9 | 49.5     | 125.7    | 55 1    | 139.7  | 45.5     | 115.6  | 45.5 11                             | 115.6  | 48 12    | 121.9   | 48 121.9   | 1.9 48   | 121.9   | 9 48  | 3 121.9 | 9 41  | 104.1  |
| Overall Width:                                   | 31.25  | 79.4  | 32.5    | 82.6  | 31.5    | 80.0  | 33 8     | 83.8     | 38 8    | 96.5   | 33 8     | 83.8   | 24.5 6                              | 62.2   | 28 7     | 71.1    | 28 71.1    | .1 28    | 71.1    | 32    | 2 81.3  | 3 25  | 63.5   |
| Overall Width<br>w/Side Table(s) Up:             | 51.5   | 130.8 | 52.5    | 133.4 | 52      | 132.1 | 53 1:    | 134.6    | 48 12   | 121.9  | 45 1     | 114.3  |                                     |        |          |         |            |          |         |       |         |       |        |
| Overall Depth:                                   | 37.5   | 95.3  | 88      | 96.5  | 33      | 99.1  | 1 44     | 111.8    | 59 1    | 149.9  | 36.5     | 92.7   | 41 10                               | 104.1  | 98       | 91.4    | 36 91.4    | .4 36    | 91.4    | 37    | 7 94.0  | 26    | 0.99   |
| Width Between Arms:                              | 21.75  | 55.2  | 22      | 55.9  | 21.75   | 55.2  | 23 5     | 58.4 2   | 25.25 6 | 64.1   | 24 (     | 61.0   | 20 5                                | 50.8   | 19 4     | 48.3    | 19.5 49.5  | .5 19.5  | 5 49.5  | 24    | 4 61.0  | 18    | 45.7   |
| Arm Height From Floor<br>(Arm Up on Drop Arm):   | 28.5   | 72.4  | 28      | 71.1  | 27.25   | 69.2  | 28 7     | 71.1     | 31 7    | 78.7   | 26.75    | 6.79   | 26 (                                | 99     | 26 (     | 92   56 | 26.75 67.9 | .9 26.75 | 6.79    | 28    | 3 71.1  | 31.5  | 08     |
| Arm Height From Floor<br>(Arm Down on Drop Arm): |        |       | 21.5    | 54.6  |         |       |          |          |         |        |          |        |                                     |        |          | 23      | 23.5 59.7  | 22.      | 5 57.2  |       |         |       |        |
| Arm Height From Seat<br>Front Section:           |        |       |         |       |         |       |          |          | 5 1     | 12.7   | 9        | 15.2   |                                     |        |          |         |            |          |         |       |         |       |        |
| Arm Height From Seat<br>Mid Section:             | 8      | 20.3  | 9       | 15.2  | 6.5     | 16.5  | 7 1      | 17.8     | 8 2     | 20.3   | 7.75     | 19.7   | 8 2                                 | 20.3   | 8 2      | 20.3    | 8 20.3     | .3 8     | 20.3    | 8.5   | 5 21.6  | 5.5   | 14.0   |
| Arm Height From Seat<br>Rear Section:            |        |       |         |       |         |       |          |          |         |        | 8.5      | 21.6   |                                     |        |          |         |            |          |         |       |         |       |        |
| Seat Height:                                     | 22.25  | 56.5  | 22.5    | 57.2  | 20.5    | 52.1  | 22   5   | 55.9     | 25.25   | 64.1   | 21 (     | 53.3   | 20 5                                | 50.8   | 20 5     | 50.8    | 22 55      | 55.9 22  | 55.9    | 20    | 50.8    | 3 27  | 9.89   |
| Seat Width:                                      | 20.5   | 52.1  | 20.5    | 52.1  | 21.5    | 54.6  | 23.5     | 59.7     | 25.75 6 | 65.4   | 23 (     | 58.4   | 19 4                                | 48.3   | 19 4     | 18.3    | 19 48.3    | .3 19    | 48.3    | 24    | 4 61.0  | 18    | 45.7   |
| Seat Depth:                                      | 19     | 48.3  | 19      | 48.3  | 19.5    | 49.5  | 22   5   | 6.53     | 24 6    | 61.0   | 20.5     | 52.1   | 20 5                                | 50.8   | 19.5 4   | 49.5    | 19.5 49.5  | .5 19.5  | 5 49.5  | 19.5  | .5 49.5 | 16    | 40.6   |
| Armrest Width:                                   | 4.5    | 11.4  | 4.5     | 11.4  | 4.5     | 11.4  | 4.5      | 11.4     | 5 1     | 12.7   | 0        | 0.0    | 2 5                                 | 5.1    | 1.5      | 11.4    | 4.5 11.4   | .4 4.5   | 11.4    | 4.5   | 5 11.4  | 2.75  | 7.0    |
| Depth Fully Reclined:                            | 71.5   | 181.6 | 71.5    | 181.6 | 65      | 165.1 | 71 1     | 180.3    | 77 19   | 195.6  | 64 1     | 162.6  | 65 16                               | 165.1  | J        | 0.0     |            | 99       | 167.6   | 9     | 0.0     |       |        |
| Depth Trendelenburg Position:                    | 75     | 190.5 | 75      | 190.5 | 75      | 190.5 | 80 2     | 203.2    | 77 18   | 195.6  | 74 1     | 188.0  |                                     |        |          |         |            |          |         |       |         |       |        |
| Back Angle<br>Upright:                           | 120°   | 120°  | 120°    | 120°  | 120°    | 120°  | 120° 1   | 120°     | 110° 1  | 110° 1 | 120°     | 120°   | 105° 10                             | 105° 1 | 110° 1.  | 110° 1  | 110° 110°  | 0. 110.  | 110°    | 110°  | 0. 110. | ٥     |        |
| Back Angle<br>With Legrest Extended:             | 130°   | 130°  | 130°    | 130°  | 130°    | 130°  | 130°   1 | 130°     | VAR \   | VAR 1  | 139°     | 139°   | 105° 10                             | 105° V | VAR V.   | VAR     |            | VAR      | NAR     | VAR   | R VAR   | ~     |        |
| Back Angle- Fully Reclined:                      | 150°   | 150°  | 150°    | 150°  | 150°    | 150°  | 150°   1 | ا50。     | VAR V   | VAR 1  | 180°   1 | 180°   | 110° 1                              | 110° 1 | 150°   1 | 150°    |            | 150°     | ° 150°  | . VAR | R VAR   | ~     |        |
| Back Height Above Seat:                          | 30     | 76.2  | 30      | 76.2  | 29      | 73.7  | 31.5     | 80.0     | 35 8    | 88.9   | 31 7     | 78.7   | 29 73                               | 73.7   | 31 7     | 78.7    | 31 78.7    | .7 31    | 78.7    | 30.5  | .5 77.5 | 16.5  | 41.9   |
| Caster Diameter:                                 | 5      | 12.7  | 2       | 12.7  | 2       | 12.7  | 5        | 12.7     | 5       | 12.7   | 4        | 10.2   | 4                                   | 10.2   | 5        | 12.7    | 5 12.7     | .7 5     | 12.7    | 5     | 12.7    | ,     |        |
| Weight:  | 128    | 58.1  | 131     | 59.4  | 131     | 58.1  | 134 6    | 8.09     | 285 13  | 129.3  | 200      | 2.06   | 80 3                                | 36.3   | 69       | 31.3    | 70 31.8    | .8 70    | 31.8    | 74    | 4 33.6  | 3 40  | 18.1   |
| Shipping Weight:                                 | 139    | 63.0  | 142     | 64.4  | 139     | 63.0  | 145 6    | 65.8     | 308 13  | 139.7  | 240 1    | 108.9  | 96 4:                               | 43.5   | 76 3-    | 34.5    | 78 35.4    | .4 78    | 34.5    | 80    | 36.3    | 3 45  | 20.4   |
| Maximium Weight Capacity<br>Evenly Distributed:  | 350    | 158.8 | 350     | 158.8 | 290     | 131.5 | 450 21   | 204.1    | 200 3   | 317.5  | 350 1    | 158.8  | 250 11                              | 113.4  | 250 11   | 113.4   | 250 113.4  | 3.4 250  | 113.4   | 4 400 | 181.4   | 350   | 158.8  |

|                      | Bariatrics               |                                  |                                      | •  |                                 |  |  | •   |                                   | •                                    |
|----------------------|--------------------------|----------------------------------|--------------------------------------|--|---------------------------------|--|--|---|-----------------------------------|--------------------------------------|
|                      | Physical<br>Therapy      | •                                | •                                    |  |                                 |  |  |   |                                   |                                      |
|                      | Rehab                    | •                                | •                                    | •  | •                               | •                                      | •  | •   | •                                 |                                      |
|                      | Sub Acute<br>Care        |                                  |                                      |  |                                 |  |  |   |                                   | •                                    |
|                      | Patient<br>Room          | •                                | •                                    | •  | •                               | •                                      | •  | •   | •                                 | •                                    |
|                      | Emergency<br>Room        |                                  |                                      |  | •                               | •                                      | •  | •   |                                   |                                      |
|                      | PACU                     |                                  |                                      | •  | •                               | •                                      | •  | •   | •                                 | •                                    |
|                      | Same Day<br>Surgery/ ASC |                                  | •                                    | •  | •                               | •                                      | •  | •   | •                                 | •                                    |
| SE                   | Post Opera-<br>tive Care |                                  | •                                    | •  | •                               | •                                      | •  | •   | •                                 | •                                    |
| SOFU                 | IV<br>Therapy            |                                  |                                      | •  | •                               | •                                      | •  | •   | •                                 | •                                    |
| PRIMARY AREAS OF USE | Pain<br>Management       |                                  |                                      | •  | •                               | •                                      | •  | •   | •                                 | •                                    |
| PRIM                 | Oncology                 |                                  |                                      | •  | •                               | •                                      | •  | •   | •                                 | •                                    |
|                      | Nocturnal<br>Dialysis    |                                  |                                      |  |                                 |  |  |   |                                   | •                                    |
|                      | Dialysis                 |                                  |                                      |  | •                               | •                                      | •  | •   | •                                 | •                                    |
|                      | РНОТО                    |                                  |                                      |  |                                 |  |  |   |                                   |                                      |
|                      | DESCRIPTION              | Lumex Three<br>Position Recliner | Lumex Hip Chair with Adj<br>Footrest | Lumex Deluxe Preferred Care<br>Recliner-Wide | Lumex Clinical<br>Care Recliner | Lumex Deluxe Clinical<br>Care Recliner | Lumex Deluxe Clinical Care<br>Recliner with Drop Arm | Lumex Deluxe Clinical Care<br>Recliner-Wide | Lumex Ortho-Biotic II<br>Recliner | Lumex Preferred Care<br>Bed-Recliner |
|                      | MODEL#                   | 574G/FR574G                      | GF4405                               | FRS6SWG                                      | FR577RG                         | FR566G                                 | FRS66DG  | FR587W                                      | FR597G/FR597P                     | AR981642                             |



Information contained herein is subject to change. The most current and complete product information can be found on our website.







© 2016, GF Health Products, Inc. All Rights Reserved.

E&J, Everest & Jennings, Graham-Field, Hausted, and Lumex are trademarks of GF Health Products, Inc.
GF Health Products, Inc. is an ISO 13485 Certified Company.

"SIMMONS Clinical Care products manufactured by GF Health Products, Inc. under license from the trademark owner, Dreamwell, Ltd., a subsidiary of Simmons Bedding Company, Atlanta, Georgia. SIMMONS is a trademark of Dreamwell, Ltd. a subsidiary of Simmons Bedding Company.

© 2016 Dreamwell, Ltd. All rights reserved.

GF1600017RevA16