



EVERYONE
KNOWS

WE WORK BEST
UNDER
PRESSURE.

+

WE CAN HANDLE
MOISTURE TOO.

ROHO

S/A

SELECTAIR®

The world's first low air loss mattress line to carry the ROHO® name.




The ROHO name has long been synonymous with the highest quality DRY FLOATATION® products. When we decided to add a low air loss mattress system to our product group, we searched for a brand that would meet the high quality standards our clients and clinicians expect from a ROHO product. *SelectAir*® is that brand, and we now manufacture the *SelectAir* line in our ROHO facility in Belleville, Illinois.

ROHO can now offer moisture control for added patient comfort and safety.

ROHO DRY FLOATATION products are known for exceptional protection against pressure, shear and friction. *SelectAir* is a product line that also guards against moisture. It is a cost-effective alternative to fully-integrated bed systems, and is ideal for clients who need special protection from the damaging effects of excess moisture on the skin associated with perspiration.

SelectAir controls moisture with wicking action. - Air gently blows underneath a vapor-permeable cover, wicking moisture away and keeping the skin at a healthy moisture level. The wicking method has several benefits over direct air flow alternatives. It is much less likely to over dry the skin (a clear contra-indication to moist wound healing). It also reduces the risk of airborne pathogens, and provides additional comfort for patients who tend to feel cold or simply find it unpleasant to have air blowing directly on them.

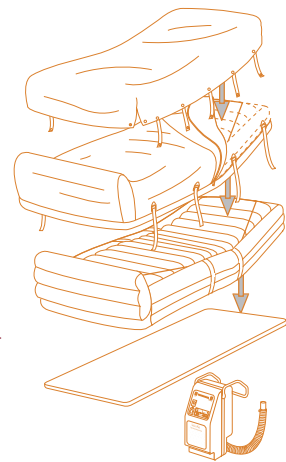

All *SelectAir* mattresses can be used on patients with any stage ischemic (pressure) ulcer, myocutaneous skin flaps, skin grafts and/or patients with temperature, moisture control, or pain management issues.




The HybridSelect™ is a low air loss mattress overlay system that revolutionizes the overlay market by providing customized comfort and optimal skin protection. HybridSelect with HealFloat™ combines the power of *SelectAir* low air loss systems with the skin protection of a ROHO HIGH PROFILE® Cushion in the heel area to create an optimal environment for skin healing and protection.

SelectAir Wicks Away Moisture Problems.

Wicking Action

Air escapes through tiny holes on the top of the mattress, and circulates under a vapor-permeable cover, to wick moisture away from patient.

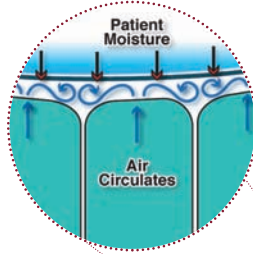

Keeps Patient Cool and Dry

"True" Low Air Loss technology allows air to escape from the top of the mattress and directly interact with moisture at the patient/bed boundary.

Exceptional Foot Protection

Foot pillow to address patient migration.



Integrated Side Bolsters

Keeps patient centered in the bed and reduces the risk of entrapment.

User-Friendly Control Unit

Rapid inflation, easy set-up, and a motor so quiet, the patient will hardly know it's there!

Quick-Change Filter

Filter is easy to access, and quick to change for hassle-free maintenance.

HybridSelect with HealFloat™ Feature Provides Ultimate Skin Protection.

HealFloat™ Feature


Provides optimal heel protection while the patient is lying in the bed.



ROHO® HIGH PROFILE® Cushion

Simply remove the ROHO HIGH PROFILE Cushion and place on any seating surface for continuous skin protection.

The ROHO Solution to Pressure, Shear, Friction, and now Moisture!

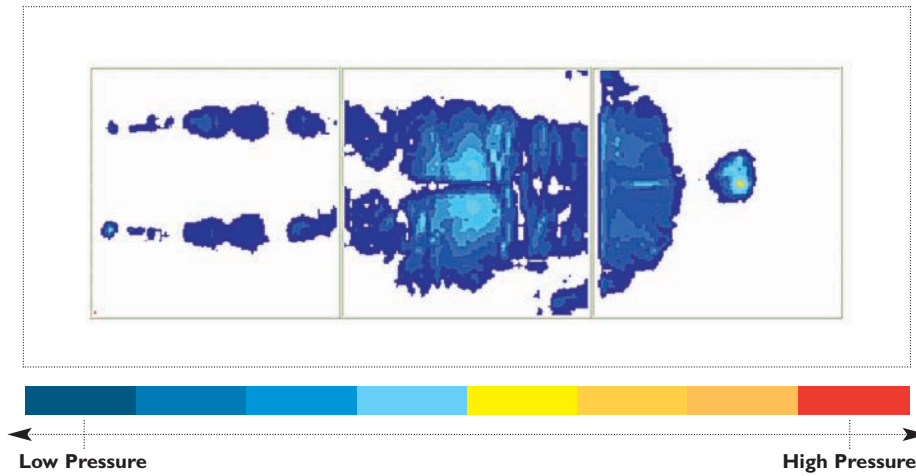

Pressure Imaging test shows the effectiveness of the *SelectAir* Mattress. Pressure is evenly distributed, leaving no areas vulnerable to ischemic (pressure) ulcers.

“True” Low Air Loss Technology


The *SelectAir* system quietly circulates air through individual cells and allows air to escape from the mattress.

Upright Patient Mode

Patient can sit upright with no loss of therapeutic benefits.



Wicking Action

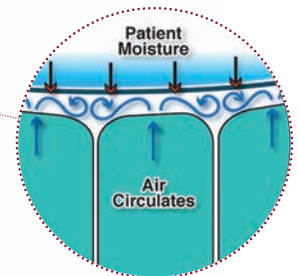

Supreme Patient Comfort

Individual cells plus side bolsters and foot pillow.

Wicking Action

Air escapes from tiny air openings, distributed equally across the top mattress surface, and circulates under a vapor-permeable cover, to wick moisture away from the patient.

SelectAir®



Quiet, portable, lightweight and very reliable, the *SelectAir* mattress and control unit provide a cost-effective alternative to other pressure management systems.



- Firm Mattress Mode: Helps facilitate nursing care during dressing changes and patient repositioning.
- Upright Patient Mode: Adds support when the head of the bed is elevated.
- Foot Pillow: Helps to prevent unnecessary skin breakdown due to patient migration.

SelectAir® Max



SelectAir Max provides all the benefits of SelectAir plus these special features:

- System Lock: Helps prevent accidental adjustments, and automatically resets the system in case of power failure.
- Auto Comfort: Takes the guesswork out of programming. During set-up, mattress sets optimal pressure according to patient's weight distribution.
- Pulsation (alternating pressure): Provides three cycle settings that continually readjust the system to ensure patient comfort.

Bariatric patients enjoy ultimate comfort on the BariSelect Mattress. It has added depth and a 680-pound weight capacity. It offers more comfort to patients, is less expensive than other fully-integrated systems, and is easier for healthcare workers to operate.



BariSelect™



- Quickset: Allows you to input the patient's actual weight, so you can preset the pressure.
- Optional Hand-Held Control: Allows patient to make adjustments as necessary.
- Added mattress depth, extended widths and specialized software accommodate extra weight capacity.

HybridSelect™



HybridSelect™ provides the benefits of SelectAir in a convenient Mattress Overlay featuring HealFloat™:

- System Lock: Helps to prevent accidental adjustments, and automatically resets the system in case of power failure.
- Firm Mattress Mode: Helps facilitate nursing care during dressing changes and patient repositioning.
- HealFloat: ROHO HIGH PROFILE Cushion inserted in the heel area creates an optimal environment for skin healing and protection.
- Pulsation (alternating pressure): Provides two cycle settings that readjust the system to ensure patient comfort.

	SelectAir	SelectAir Max	BariSelect	HybridSelect
Mattress Configurations	Individual baffles divided into 3 sections; head, seat and foot. Integrated side air bolsters with foot pillow.	Individual baffles divided into 3 sections; head, seat and foot. Integrated side air bolsters with foot pillow.	Individual baffles divided into 3 sections; head, seat and foot. Integrated side air bolsters with foot pillow.	Individual baffles divided into 3 sections; head, seat and foot. Integrated side air bolsters. ROHO HIGH PROFILE cushion in foot section.
Weight	Control Unit – 9 lbs. (4.0 kg) Mattress – 10 lbs. (4.5 kg)	Control Unit – 9 lbs. (4.0 kg) Mattress – 10 lbs. (4.5 kg)	Control Unit – 9 lbs. (4.0 kg) Mattress – 15 lbs. (6.8 kg)*	Control Unit – 9 lbs. (4.0 kg) Mattress Overlay – 5 lbs. (2.3 kg) Cushion – 3.5 lbs (1.6 kg)
Voltage	120 V AC 50/60 Hz	120 V AC 50/60 Hz	120 V AC 50/60 Hz	120 V AC 50/60 Hz
Current	Less than 5 amps	Less than 5 amps	Less than 5 amps	Less than 5 amps
Certification	UL	UL	UL	UL
Mattress Weight Limit	400 lbs. (181.6 kg)	400 lbs. (181.6 kg)	680 lbs. (308.7 kg)	300 lbs. (136 kg)
Mattress Sizes Available	Fits most standard bed frames. Custom sizes available.	Fits most standard bed frames. Custom sizes available.	Fits most bariatric bed frames. Custom sizes available.	Fits most standard bed frames.
Cover (Top Surface)	Waterproof, vapor permeable, flame resistant and biocompatible zippered mattress cover	Waterproof, vapor permeable, flame resistant and biocompatible zippered mattress cover	Waterproof, vapor permeable, flame resistant and biocompatible zippered mattress cover	Waterproof, vapor permeable sheet
Warranty	One Year Limited Warranty on control unit and soft goods	One Year Limited Warranty on control unit and soft goods	One Year Limited Warranty on control unit and soft goods	One Year Limited Warranty on control unit and soft goods; Two-Year Limited Warranty on removable HIGH PROFILE cushion.

* Based on 48" mattress

Note: U.S. Medicare Part B Code E0277 and E0372 (HybridSelect). Please note that decisions on how to code a particular product are solely the responsibility of the Medicare Provider (i.e., the supplier). In addition, Medicare coding and coverage criteria are subject to change. As such, The ROHO Group strongly recommends that a supplier should contact their respective DMERC and/or SADMERC to verify current coding and coverage criteria needed to qualify for reimbursement under a specific code.



LITG00623, SCHW, 5M, 04/07
©2007 The ROHO Group.

ROHO®, DRY FLOATATION®, SelectAir®, and BariSelect™, HybridSelect™, HealFloat™, HIGH PROFILE® Shape-Fitting Technology® and "The ROHO Group Design™" are registered trademarks of ROHO, Inc. The ROHO Group, Inc. has a policy of continual product improvement and reserves the right to amend specifications presented in this brochure.